

# Platinum 6000™

Low Air Loss Alternation/Immersion Therapy



Shown with Rev. 2 control unit  
and navy top cover

Platinum 6000™ is ideal for facilities wanting ease of operation and effective microclimate management. Cell pressure is measured twice per second to maintain optimal patient immersion and envelopment. Preprogrammed comfort and therapy settings are based on patient height and weight. Can be sized to fit most bed frames.

- *Two therapy modes (alternation and immersion) for pressure injury prevention and treatment*
- *First Point of Contact™ top cover free of harmful chemicals*
- *Transport mode keeps mattress inflated without power*
- *Effective in managing patient microclimate*


sizewise.com 800-814-9389

  
**Sizewise**®  
Designed to heal

1



2

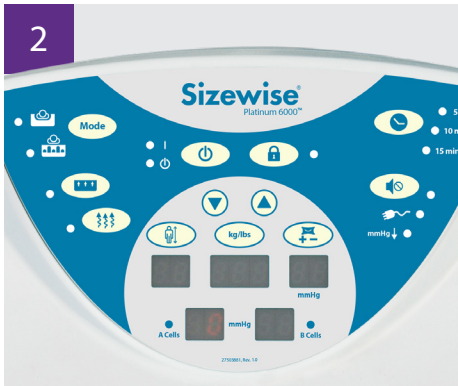

3



4



- 1 — Alternation and immersion therapies for prevention and treatment of pressure injuries
- 2 — Preprogrammed comfort and support settings based on patient height and weight (Rev. 2 keypad shown)
- 3 — Transport mode keeps mattress inflated without power
- 4 — Available for a variety of facility-owned or Sizewise rental bed frames

## Key Specs

### MATTRESS

Capacity ..... 600 lb./900 lb.  
 Width ..... 35", 39", 42", or 48"  
 Length ..... 80" or 86"  
 Standard and Bariatric Sizes ..... Available  
 Mattress Overlay ..... Available

### CONTROL UNIT

Liters of Air Per Minute ..... 25  
 Zoned Low Air Loss ..... Yes  
 Alternation Cycles ..... 5, 10, or 15 Min.  
 Transport Mode ..... Yes  
 Alarms ..... Power Failure and Low Pressure  
 Pressure Sensor Technology ..... Yes  
 Unit Weight ..... 12 lb.

### FABRIC OPTIONS

Extreme Polycarb 4-Way Stretch ..... Standard  
 Non-Skid Polyurethane Bottom Cover ..... Standard

### FABRIC CHARACTERISTICS



Powered Surface



4-Way Stretch



Durability

sizewise.com 800-814-9389

