


Pediatric Pulse™

Low Air Loss Pulsation Therapy

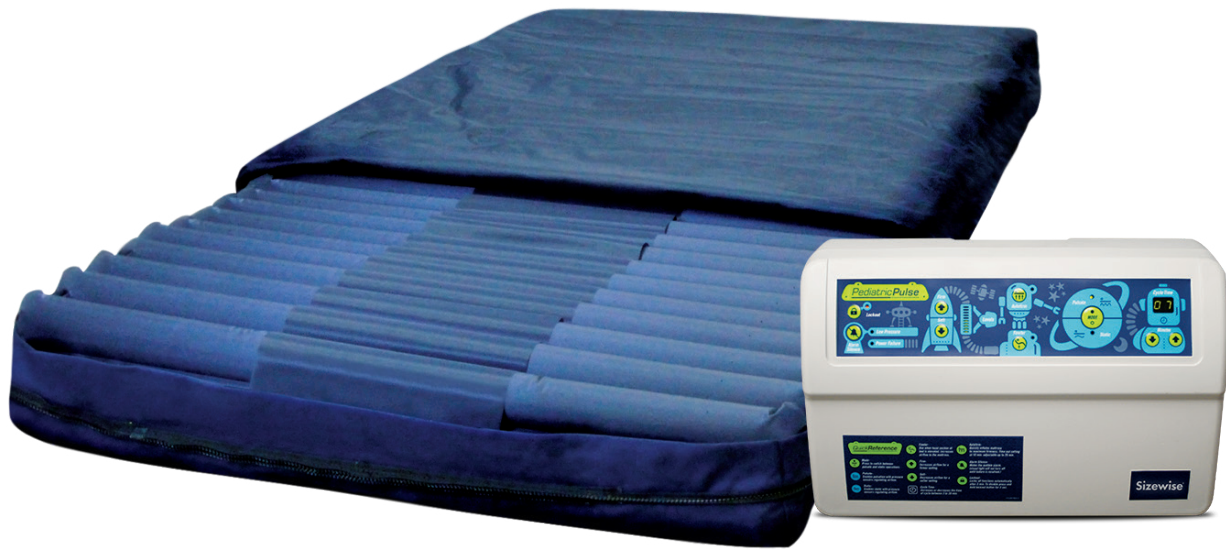


Pediatric Pulse™ is a latex-free mattress replacement system designed for pressure injury treatment and prevention in children. The surface's soothing pulsation offers a nonpharmaceutical treatment for Neonatal Abstinence Syndrome (NAS). The top cover is vapor-permeable, low-friction, and anti-shear. The three-zoned surface features vented air cells to provide pulsation or static therapy and manage microclimate. Radiolucent materials for ease of X-ray.

- Adjustable pulsation times from three to 20 minutes
- Vapor-permeable, low-friction, anti-shear top cover
- Kid-friendly control unit keypad design
- Autofirm feature for safe and comfortable transport
- Available in standard and custom sizes

size-wise.com 800-814-9389


Size-wise®
Designed to heal

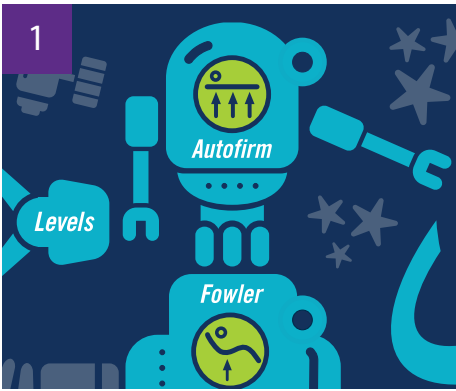


Key Specs

MATTRESS

Minimum Patient Weight 1.6 lbs.
 Maximum Patient Weight 100 lbs.
 Sizes Standard and Custom

CONTROL UNIT

Alarms Power Failure and Low Pressure
 CPR Release Yes
 Pressure Sensor Technology Yes
 Pulsating Low Air Loss 3-20 Min.
 Unit Weight 11 lbs.
 Zoned Low Air Loss Yes

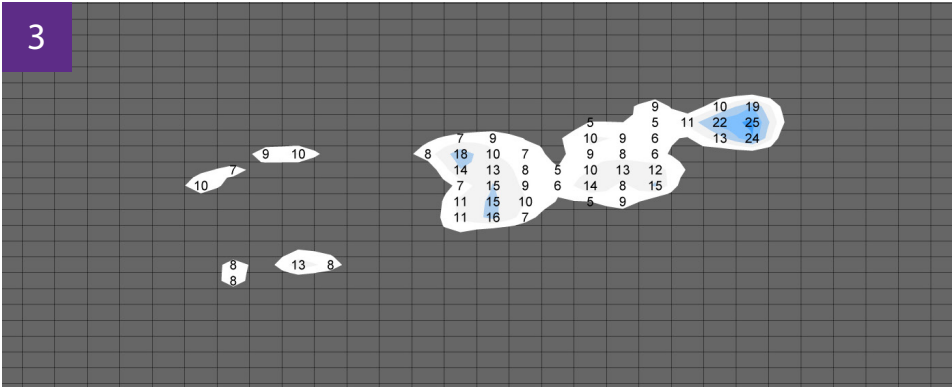

- 1 — Autofirm feature for safe and comfortable transport
- 2 — Easy-to-use controls and quick-glance instructions
- 3 — Effective PI prevention and treatment for patients as small as 1.6 lbs. (pressure mapping features male patient, 3 yo, 34.8 lbs.)

sizewise.com 800-814-9389

