


Carewise 7.1™

Foam Support Surface

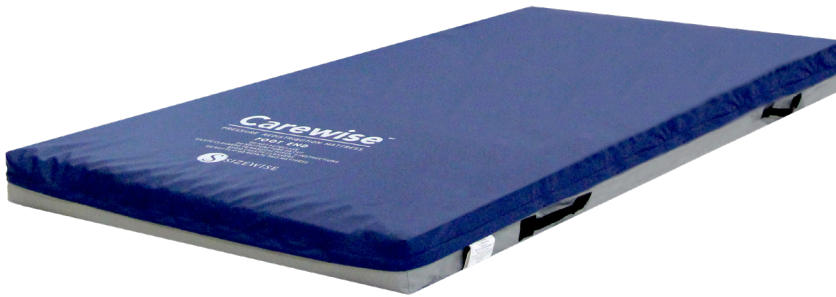

Key Specs

Capacity	500 lbs.
Width	35", 39", or 42"
Length	80" or 84"
Height.....	7"

CONSTRUCTION

Foam Layers.....	4
Defined Perimeter.....	Optional

FABRIC OPTIONS

4-Way Stretch Top Cover	Standard
PVC Laminate Bottom Cover	Standard
Fire Barrier	Standard

Carewise 7.1™ features a viscoelastic foam layer for pressure redistribution, shear reduction, and air circulation, keeping the surface cool for optimal patient comfort. The heel area is softer than the rest of the surface to protect vulnerable heels, and a concealed zipper prevents skin tears and protects inner foam from fluid ingress.

- *Diamond-cut top layer to keep patient cool*
- *Heel area is softer than rest of surface*
- *Concealed zipper to prevent skin tears and protect inner foam*
- *Two kinds of viscoelastic foam to "float" patient on common pressure points*

size-wise.com 800-814-9389

