

# Envy Line® O Series™

OR Table Pad




Envy Line O Series features F3 Free™ construction that reduces harmful chemicals in the fabric, foam, and fire barrier. A four-way stretch top cover allows immersion into the pad for optimal pressure redistribution in order to prevent pressure injuries. The universal design is radiolucent and fits most OR tables.

- *F3 Free construction reduces harmful chemicals in the fabric, foam, and fire barrier*
- *Optimal immersion for pressure redistribution*
- *Five-piece universal design fits most OR tables*
- *Custom sizes available*

sizewise.com 800-814-9389

  
**Sizewise®**  
Designed to heal

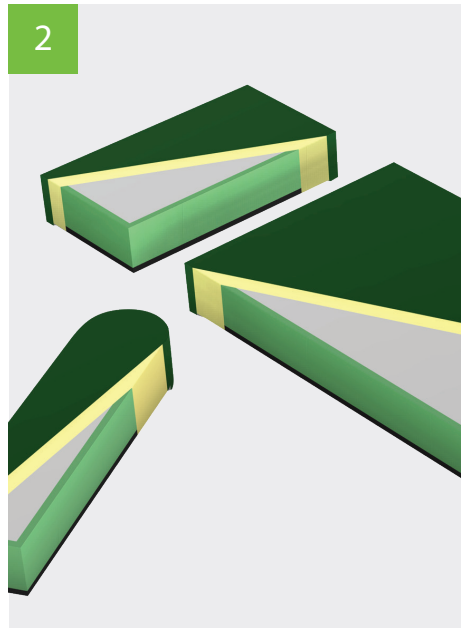

## Key Specs

### BED

|                    |                 |
|--------------------|-----------------|
| Capacity .....     | 500 lbs.        |
| Head .....         | 20"W x 10.5"L   |
| Torso .....        | 20"W x 42.375"L |
| Leg .....          | 20"W x 18.5"L   |
| Arm .....          | 6"W x 26"L      |
| Custom Sizes ..... | Available       |

### FABRIC OPTIONS

|  |          |
|--|----------|
| Polycarbonate 4-Way Stretch Top Cover . . .  | Standard |
| Dual-Sided Polyurethane Bottom Cover . . . . | Standard |
| Halogen-Free Fire Barrier . . . . .          | Standard |
| MAP Liner™ . . . . .                         | Optional |

- 1 — Polymer gel topper contours to patient body while providing low-interface pressure
- 2 — EcoEase™ pieces are easily replaceable by section



A PRACTICE GREENHEALTH PROGRAM

Sizewise proudly supports Practice Greenhealth's Healthier Hospitals Safer Chemicals Challenge for Healthy Interiors. Our Envy Line products meet or exceed Healthier Hospitals goals of reducing patient and caregiver exposure to chemicals of concern.

[sizewise.com](http://sizewise.com) 800-814-9389



**Sizewise**®  
Designed to heal